

Puritan® Fecal Opti-Swab® Collection & Transport System.

Optimize your specimen collection and preservation of fecal enteric bacteria with Puritan® Fecal Opti-Swab®. A liquid base ready-to-use system for rectal swabbing or direct collection from fecal sample.




HOW IT WORKS

Simple to use – eliminates unnecessary specimen handling steps.

- Our HydraFlock® swab assures efficient collection and elution of sample
- Only a small amount of sample is needed, making handling and transporting easy
- Versatile system allows for use with traditional plating methods or liquid based automation platforms

TESTING APPLICATIONS

While different types of bacteria can cause enteric infections routine stool cultures can generally be used to screen for public health issues such as *Salmonella spp.*, *Shigella spp.*, and *Campylobacter spp.* Can also be used to culture: *E.coli*, *V. parahaemolyticus*, *Escherichia coli* O157:H, *Y. enterocolitica*, *Shigella flexneri*, *C.jejuni*, *S.typhimurium*, *C.difficile*, and *Efaecalis (VRE)*.

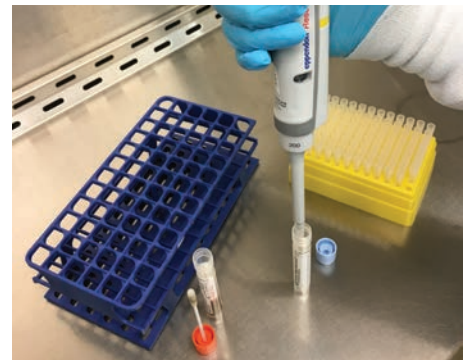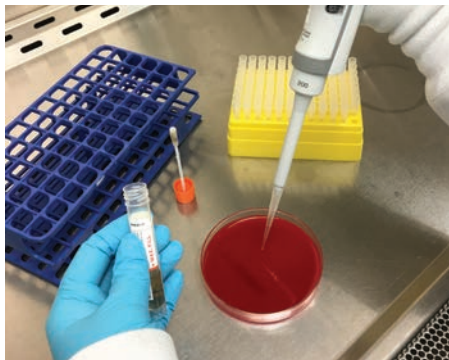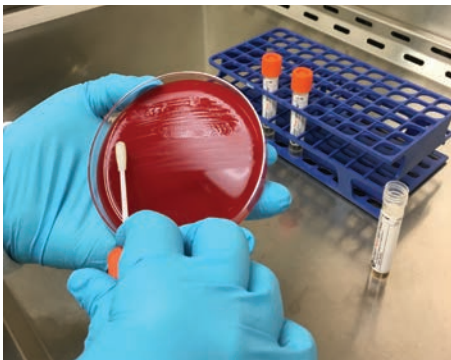

Optimize your specimen collection.



CONVENIENT, ALL-IN-ONE KIT.

Each System is enclosed in a sterile, easy-peel pouch.

LEAK-PROOF, SELF-CENTERING CAP

Self-centering cap is designed to facilitate easy swab removal and disposal.

LIQUID MEDIUM

2 ml Cary-Blair medium is compatible with manual and automated specimen processors.

MEDICAL GRADE MATERIALS

The vial allows for safe collection and transport of specimens.

MAX FILL

Indicator line showing maximum sample quantity.



HIGH ABSORPTION
Puritan® HydraFlock® provides high absorption.

MOLDED BREAKPOINT
Handle is sturdy and flexible. Offers a molded breakpoint with printed indicator line for safe insertion.



PURITAN® EXCLUSIVE
HydraFlock® is a patented proprietary manufacturing design process offered only by Puritan®.

This creates a three-dimensional microstructure fiber that enhances collection and release capabilities.

PURITAN® FECAL OPTI-SWAB® COLLECTION & TRANSPORT SYSTEM

Puritan Item No.	Product Description	Volume	Tube length x diameter	Packaging pack, weight
CB-206	Media: Cary-Blair medium Swab: One flock-tipped collection swab (HydraFlock® 3506-H)	2 ml	12 x 80 mm	6/50, 9 lbs

FOR SHIPPING AND LOGISTICS INQUIRIES, VISIT:

beckmancoulter.com/contact

FOR PRODUCT INQUIRIES, CONTACT:

+1 207-876-3130 (toll free) | sales@puritanmedproducts.com

FOR MORE INFORMATION, VISIT BECKMANCOULTER.COM/PURITAN

HydraFlock® is a registered trademark of Puritan Medical Products Co. LLC. ISO 13485:2016

Not all products are available in all countries.

© 2019 Beckman Coulter, Inc. All rights reserved. Beckman Coulter, the stylized logo, and the Beckman Coulter product and service marks mentioned herein are trademarks or registered trademarks of Beckman Coulter, Inc. in the United States and other countries.

For Beckman Coulter's worldwide office locations and phone numbers, please visit www.beckmancoulter.com/contact

FL-239790

